



## PF 300

**RELIABLE AND STABLE TEMPERATURE IS OBTAINED THROUGH AIR-FORCED COOLING.** This freezer offers a great overview of the stored items due to the double pane glass door.



AVAILABLE IN  
ATEX



AIR FORCED  
COOLING



FILTERLESS  
CONSTRUCTION

### STAINLESS STEEL & SPARK FREE

The unit is available with stainless steel AISI 304 Exterior, so that the interior and exterior is made entirely of stainless steel, providing a stronger more resistant unit.

This unit can be made in an ATEX version – approved according to EN 60079-15 category 3, zone 2. CE EX 11-/3 G

### FLEXIBLE STORAGE SOLUTIONS


Arctiko's PF range is a highly customizable allowing you to choose between a wide range of universal shelves and drawers.

This allows the user to organize their unit exactly as desired. Our universal shelf and drawer guides allow you to mix shelves and drawers for flexible storage.

### GREAT COOLING PERFORMANCE

Air forced cooling provides superior temperature uniformity and stability. Forcing the airflow around the shelves and the stored items, ensures that all parts of the unit is equally cold and avoids any cold or hot spots. Distributed by a powerful fan and open air flow areas along the sides, the air is always distributed evenly throughout the fully loaded unit.

### MAIN DIMENSION FOR PF 300

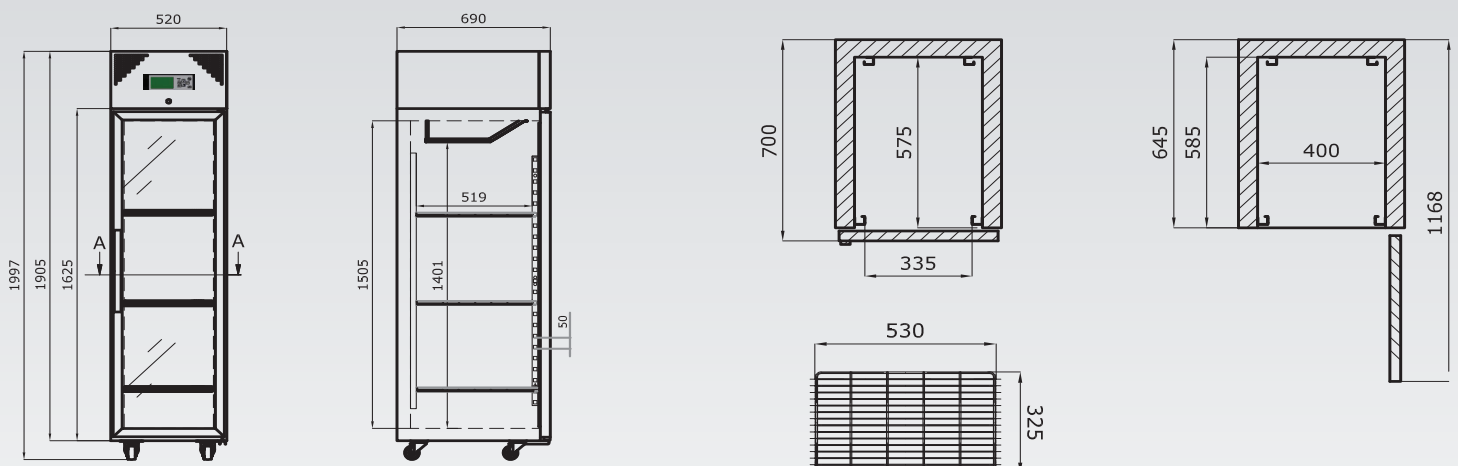

## MAIN SPECIFICATIONS

|                                |                     |
|--------------------------------|---------------------|
| Item code                      | DAI 0890 / DAI 0891 |
| Temperature range (°C)         | -10 / -23           |
| Max. ambient temperature (°C)  | 25                  |
| Capacity (L)                   | 352                 |
| Exterior dimensions WxDxH (mm) | 520x690x1997        |
| Interior dimensions WxDxH (mm) | 400x585x1505        |
| Weight (KG)                    | 103                 |
| Insulation (mm)                | 60                  |
| Cooling technology             | air forced          |
| Controller model / probe type  | G214 / PT1000       |

## POWER

|                             |       |
|-----------------------------|-------|
| Power supply (V)            | 230*  |
| Frequency (Hz)              | 50/60 |
| Power consumption (kWh/24H) | 11,9  |
| AMP (A)                     | 4,0   |
| Power (Watt)                | 613   |
| Horse power (HP)            | 1/2   |

## ALARM FUNCTIONS

|                                          |           |
|------------------------------------------|-----------|
| Alarm display as text - not codes        | yes       |
| Visual / acoustic alarm                  | yes / yes |
| Power failure alarm                      | yes       |
| Adjustable high / low temperature alarm  | yes / yes |
| Open door alarm                          | yes       |
| Probe failure alarm                      | yes       |
| Low battery alarm                        | yes       |
| Compressor failure alarm                 | no        |
| Battery backup for alarms (approx. hrs.) | 72        |

## DATA LOGGING & EXTERNAL CONNECTION

|                                                               |          |
|---------------------------------------------------------------|----------|
| Voltage free contact for remote alarm (e.g. GSM alarm module) | yes      |
| Electronic data logger                                        | yes      |
| RS 485 port                                                   | yes      |
| USB port for software update & data download                  | yes      |
| Temperature chart recorder                                    | optional |

## DISPLAY FEATURES

|                                               |                         |
|-----------------------------------------------|-------------------------|
| Display type                                  | 4,3" with touch buttons |
| Number of optional reference probes           | 2                       |
| Password protection for turning unit on / off | yes                     |
| Temperature graph on display                  | yes                     |

## OTHER FEATURES

|                                                                   |          |
|-------------------------------------------------------------------|----------|
| Automatic defrost                                                 | yes      |
| Antifreeze function                                               | no       |
| Vacuum release port                                               | no       |
| Access port for external probe etc. (int. diameter mm)            | 12       |
| Light                                                             | yes      |
| Lock with key                                                     | yes      |
| Castors (total pcs) / castors with brake (pcs)                    | 4 / 2    |
| Fixed feet (pcs) / levelling feet (pcs)                           | no       |
| Sub lids / inner doors (pcs)                                      | no       |
| Reversible door(s)                                                | yes      |
| Fixed shelves (pcs) / adjustable shelves (pcs)                    | 0 / 3    |
| Max. no. of shelves (pcs) / recommended max. no. of shelves (pcs) | 23 / 11  |
| Max. load per shelf (kg)                                          | 30       |
| Adjustable drawers (pcs)                                          | optional |
| Max. no. of adjustable drawers (pcs)                              | 8        |
| Max. load per adjustable drawer (kg)                              | 30       |

## REFRIGERATION & CABINET

|                             |                                 |
|-----------------------------|---------------------------------|
| Number of compressors       | 1                               |
| Refrigerant                 | R290                            |
| Interior cabinet material   | stainless steel                 |
| Exterior cabinet material   | painted steel / stainless steel |
| Colour (exterior cabinet)   | white / stainless steel         |
| Foam type                   | polyurethane                    |
| VIP vacuum insulated panels | no                              |

## SHIPPING

|                                |           |
|--------------------------------|-----------|
| Shipping weight (KG)           | 160       |
| Shipping dimensions WxDxH (cm) | 83x67x220 |
| Shipping volume (m³)           | 1,22      |

\*110V/60Hz available as option

## OPTIONAL ACCESSORIES

GSM alarm module  
 Temperature chart recorder  
 Shelves (perforated, wire plastic coated, wire stainless steel)  
 Drawers (stainless steel)  
 Modbus adapter  
 Voltage stabilizer  
 Voltage converters  
 LED light

## Arctiko A/S

Oddesundvej 39 · DK-6715 Esbjerg N · Phone: +45 7020 0328  
 info@arctiko.com ■ arctiko.com

