

Dog/Large Animal Run

MACS[®]

MODULAR ANIMAL CAGING SYSTEM




**Build With This
Modular Design and SAVE!**



ALTERNATIVE DESIGN
Lab Animal Housing SolutionsSM

Select the Options That Meet Your Needs

Door Options

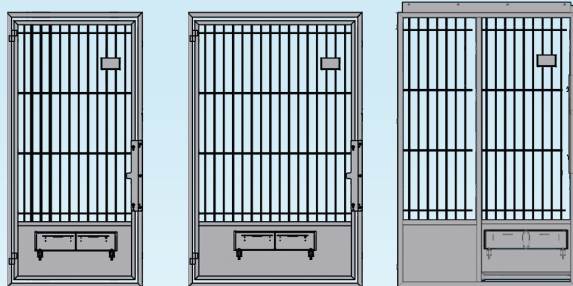


3 Foot Swing-out Door
34.75" W x 74" H

4 Foot Swing-out Door
46.75" W x 74" H

4 Foot Sliding Door
46.75" W x 76.625" H


- Standard front panels are 3' & 4' swing-out door or 4' sliding door.
- Doors are full rod in square tube frame. Standard rod spacing is 2.5".
- Sliding & swing doors have a gravity latch close system. Doors can be padlocked closed (lock not supplied).

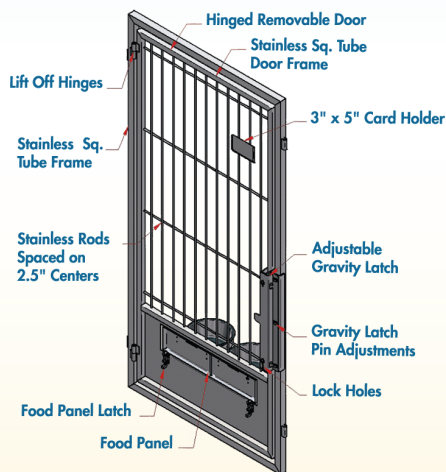

2 Qt. Food/Water Bowls (Standard)




3 Qt. Food/Water Bowls (Option)

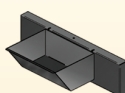

Feed Trough for Swine (Option)

- Food panel door swings around for access to feed and water bowls.
- Door accommodates 2 two quart bowls for feed and water as standard. Three quart bowls, feed trough and automated water are available options.


Side Options

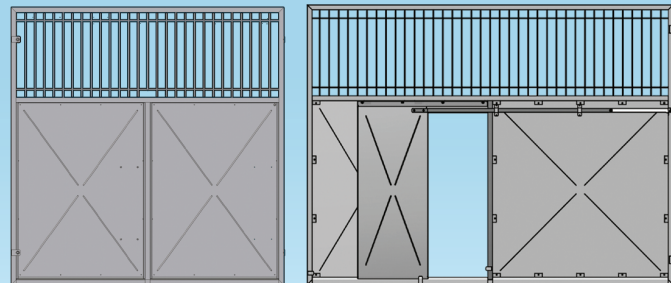

6 Foot Stainless Steel Panel
72" W x 74" H


8 Foot Stainless Steel Panel With Social Door
96" W x 74" H

- Standard side panel sizes are 6 foot and 8 foot deep.
- Side panels have stainless steel (Standard) or Trespa®* (Option) solid bottom and bar top edge for more light. Top 2' is rod construction, bottom 4' is solid. Standard rod spacing is 2.5".
- All 304 stainless steel tube frame construction. Solid SS sheet panels and Trespa panels are bolted in.
- Socialization door is standard. Solid panel is used at run ends.

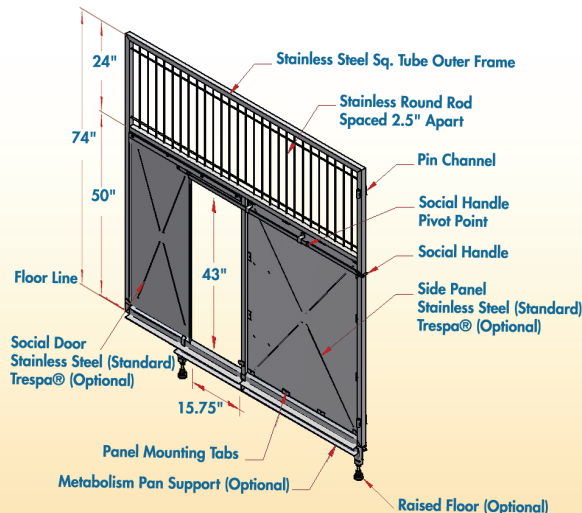

Select floor Option


- Panels that sit on facility floor are standard.
- Panels for raised floor unit.
- Leveling foot on each leg

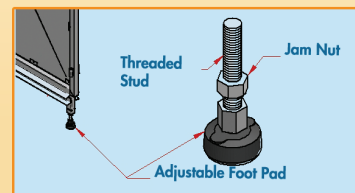


Raised Floor Option

*Standard color for Trespa® is Heaven Blue.

Back Options

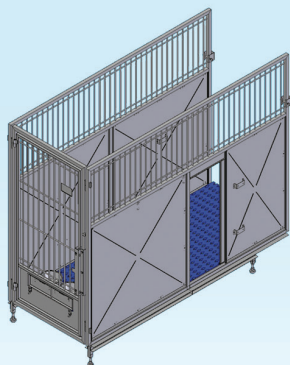


Pen Back Attached to Facility Wall (Standard)

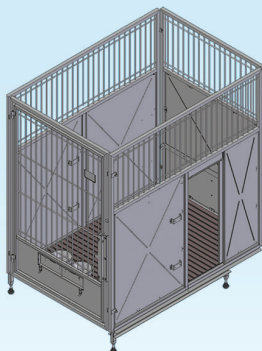

Back Panel (Option)

Pens may be free-standing
(with the addition of a back panel)

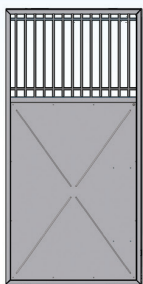


Back Panel
34.75" W x 74" H

- Back panels are 3' and 4' wide.
- Top 2' is rod construction, bottom 4' is solid SS sheet panel. (Standard)
- Trespa solid panels may be added as an option. Solid SS sheet panels and Trespa panels are bolted in.


Assembly Options

Front/side/back panels

- Pin together (Standard)
- Bolt together (Option)

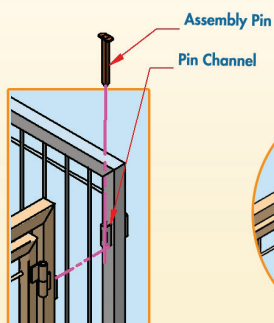


Pin Assembly (Standard)

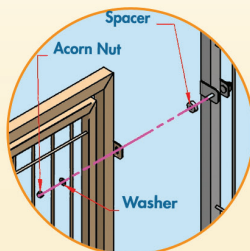

Bolt Assembly (Option)

Floor Options


Runs with raised flooring must be selected in order to have flooring options.

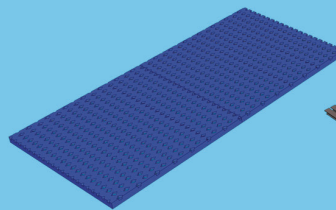

Tenderfoot® (Standard)

Tenderfoot® (Standard).

- Flooring is 1" #9 Tenderfoot. Blue plastisol coating over steel. Orange plastisol is an option.
- The run floor is supplied in sections for easier technician handling.
- This floor can support up to a total animal weight of 250 lbs.

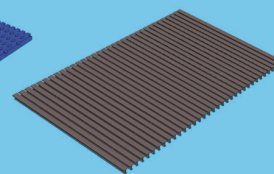


Fiberglass Slat (Option)

Fiberglass Slat (Option)


- 1-5/8" wide gray fiberglass slat with 1/2" gap.
- The run floor is supplied in sections for easier technician handling.
- This floor can support up to a total animal weight of 250 lbs.

Other options on back

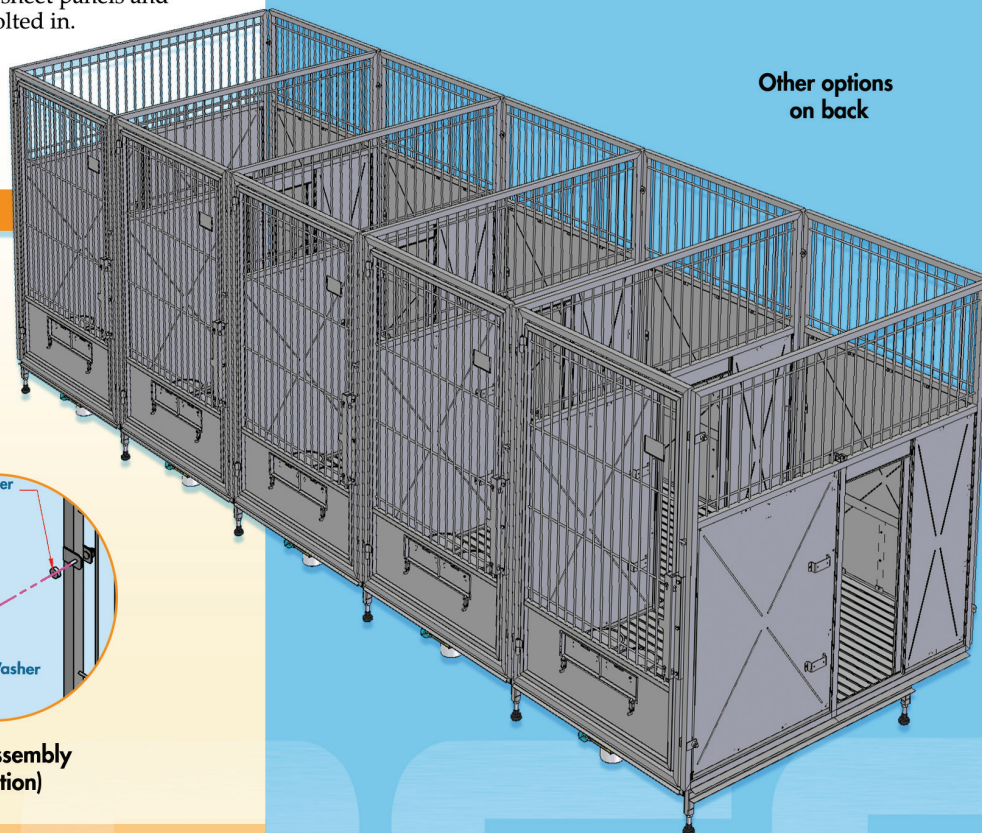


MACS
MODULAR ANIMAL CAGING SYSTEM®

Other Options

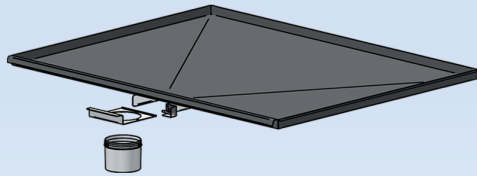


Metabolism Pan (front half of run)

- On raised floor units a metabolism pan can be added as an option. A divider panel is required to keep the dog in the front half of the pen when metabolism study is being performed.
- A 32 oz. collection jar is included which slides forward for easy access.

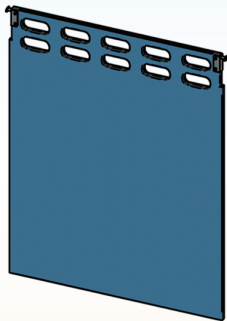


Divider Panel

- The divider panel is formed from 1/2" Trespa® sheet. Standard color is Heaven Blue.
- The panel slides into position from the top and is held in place via pins and guides on each vertical side.
- Stainless steel panels are available as an option.
- Size is 4' high by 3' or 4' wide.

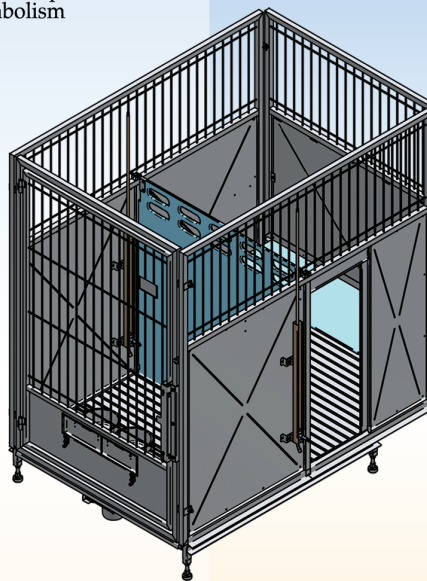

Automated Water System by Edstrom or SE

- Stainless steel Dog/Large Animal Manifold with adjustable mount bracket.
- Universal QD plug inlet.
- Kennel Divider/Panel Mount.
- Drinkers available for: dogs, hogs, sheep, goats.

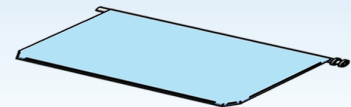

Resting Bench

- The resting bench is approximately 30.5" wide by 12.5" deep and made of Trespa.
- The bench sits 14" from the floor.
- Standard color is Heaven Blue.
- The bench folds down when not needed for more space.
- The bench attaches to the side panel with appropriate fasteners.

Warranty

Alternative Design warranties all materials* and workmanship for a period of one year. Should failure of unit occur, Alternative Design will, at its option, repair or replace the defective components. This warranty excludes damage due to abuse, misuse, as well as normal wear from use.

*Automated water is covered under the original equipment manufacturer's warranty for a period of one year.

www.altdesign.com

1.800.320.2459

Modular Animal Caging System® is a registered trademark of Alternative Design Manufacturing and Supply, Inc.

Trespa® is a registered trademark of Trespa International B.V.

Tenderfoot® is a registered trademark of Tandem Products, Inc.

ALTDESIGN

COPYRIGHT © 2007 ALTERNATIVE DESIGN
MANUFACTURING AND SUPPLY, INC.

Printed in the USA

Version A8-07