



Small & Large Animal



**Alternative Design**  
MANUFACTURING & SUPPLY

*Proven Products. Creative Solutions.*

# Small & Large Animal Products



## Housing

|                                                 |          |
|-------------------------------------------------|----------|
| <b>Canine Tandem Unit</b>                       | <b>3</b> |
| <b>Ferret Housing Unit</b>                      | <b>3</b> |
| <b>Modular Animal Runs</b>                      | <b>4</b> |
| <b>Non-Human Primate Enrichment Center Unit</b> | <b>6</b> |
| <b>Non-Human Primate Group Housing Unit</b>     | <b>6</b> |
| <b>Non-Human Primate Quad Unit</b>              | <b>7</b> |
| <b>Rabbit Housing 6-Cage 4 Sq. Ft. Unit</b>     | <b>8</b> |
| <b>Rabbit Housing 6-Cage 6 Sq. Ft. Unit</b>     | <b>8</b> |
| <b>Tree Shrew 2-Cage Unit</b>                   | <b>9</b> |

## Isolation & Containment

|                                   |          |
|-----------------------------------|----------|
| <b>Glove Port Animal Isolator</b> | <b>9</b> |
|-----------------------------------|----------|

## Transport Units

|                                         |           |
|-----------------------------------------|-----------|
| <b>Animal Holding Transport Unit</b>    | <b>10</b> |
| <b>Open Top Animal Transport Unit</b>   | <b>10</b> |
| <b>Animal Transport Cart with Scale</b> | <b>11</b> |

For more information about the features and specs, including part numbers and product dimensions on each individual product, scan the QR code to visit our website.

Overview


Features

Specs

Care

Warranty



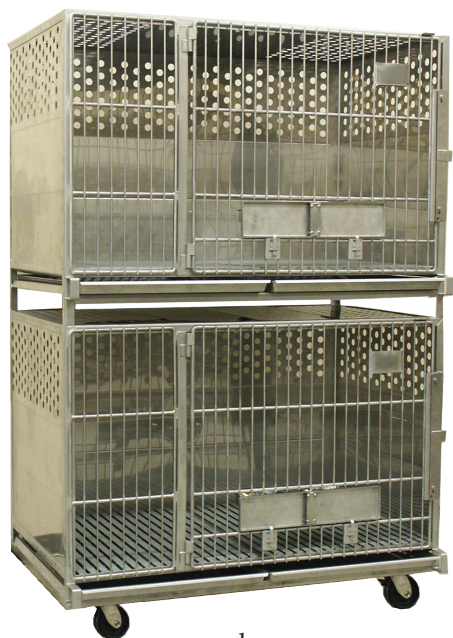

scan here



more info

## Standard Features Include:

- Frame made from 304 stainless steel
- 304 stainless steel cage sides and top with perforated holes for more light
- 5" casters, 2 with brakes
- Large cage access door with swing out feed and water bowls
- Gravity latch for positive closing with padlock capability
- Cardholder welded onto front door

## Part Numbers:

LGA2-40X34X82 (8.6 Sq. Ft. Floor)  
LGA2-55X34X82 (12.0 Sq. Ft. Floor)

## Dimensions:

40"W x 34.37"D x 82"H  
55" W x 35 3/8" D x 82" H

Designed for use with:



canine

## Ferret Housing Unit



scan here



more info

## Standard Features Include:

- Made from 304 stainless steel including the frame, feeders and pans
- Autoclavable
- Side and back panels have solid lower with wire rod upper section for more light and visibility with better ventilation
- J feeder and 32 oz. external mount water bottle
- Removable polypropylene floors for animal comfort and quieter housing
- Clean or disinfect up to 180° F

## Part #

FT6-18X20X18-WB

## Dimensions:

45" W x 24" D x 75" H

Designed for use with:



ferret



## Standard Features Include:

- Standard front panels are 3" & 4" swing-out door or 4" sliding door
- Doors are full rod in square tube frame. Standard rod spacing is 2.5"
- Socialization door is standard. Solid panel is used at run ends
- All 304 stainless steel tube frame construction. Solid SS sheet panels and Trespa panels are bolted in
- Standard side panel sizes are 6 foot and 8 foot deep
- Facility flooring is the standard installation

## Designed for use with:



canine



swine



sheep



goat



# MACS

## MODULAR ANIMAL CAGING SYSTEM

Exclusively from  
**ALTERNATIVE DESIGN**

# A Better Fit for Your Research

**MACS**<sup>®</sup>  
MODULAR ANIMAL CAGING SYSTEM



## Large Animal Run

Create modular run sizes of 18, 24 or 32 sq ft to maximize your space

Customize the run by selecting the options that best meet your needs


Simplify future change or expansion with our modular design

**5** Year  
Warranty  
Stainless Steel  
Fabricated Products

Designed to house:



Canine

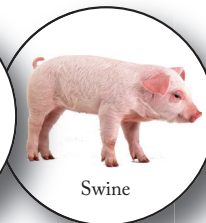


Swine



Sheep

### Options to choose from:

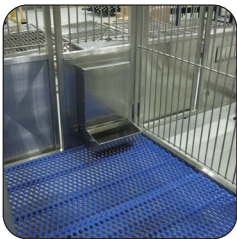

Floors



Socialization




Metabolism



Socialization



Feed/Water

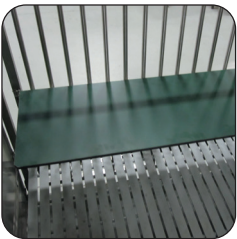


Resting Bench



Doors



Feed

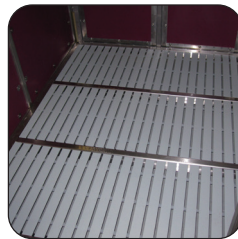

Floors



Doors

**Alternative Design**  
MANUFACTURING & SUPPLY

*Proven Products. Creative Solutions.*

[www.altdesign.com](http://www.altdesign.com)

REQUEST A QUOTE TODAY

☎ 800-320-2459

✉ [info@altdesign.com](mailto:info@altdesign.com)





## Standard Features Include:

- Made from 304 stainless steel construction
- 4-section design
- Each section has casters so sections can be separated for transport and cleaning
- Each section has 2 locking brakes
- Center section has perches and benches at various levels
- Artificial turf for foraging as well as foraging trays and feeders

**Part #**

PR6.2-128X60X82

**Dimensions:**

128" W x 60" D x 82" H

Designed for use with:



non-human primate

## Non-Human Primate Group Housing Unit



## Standard Features Include:

- Balcony – allows animals to see down hallways and aisles
- Door latch pin
- Treat door – allows placement of treats or toys into the enclosure without staff entering or disrupting the animals
- Removable panels allow for optional socialization
- Perch or bench – allows use of vertical space and a pathway to the balcony
- Holding cages – up to 4 per pen

**Dimensions:**

74.5" W x 66.5" D x 95.5" H

Designed for use with:



non-human primate



#### Standard Features Include:

- Made from 304 stainless steel construction
- External mount water bottles with automated water system option
- 2-handle easy gliding squeeze back with directional baffle
- Welded excreta pan
- J feeder
- Tamper resistant latches and card holders
- Built-in 4-level horizontal socialization

#### Part Numbers:

PR4C6.3 (6.3 sq. ft. floor)  
PR4C4.5 (4.5 sq. ft. floor)

#### Dimensions:

72" W x 36.5" D x 83" H  
65.5" W x 31.5" D x 83" H

Designed for use with:




non-human primate



# Primates are social climbers

*It's Only Natural :*

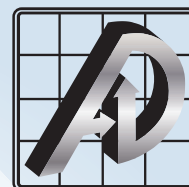

## Nonhuman Primate Group Housing

*Designed to meet the spatial and enrichment needs of primates as described in Euro ETS 123 and new ILAR Guide.*

For more details, go to:

[www.altdesign.com](http://www.altdesign.com)



scan here



more info

## Standard Features Include:

- Made from 304 stainless steel including the frame, feeders and pans
- Autoclavable
- Side and back panels have solid lower with wire rod upper section for more light and visibility with better ventilation
- J feeder and 32 oz. external mount water bottle
- Removable polypropylene floors for animal comfort and quieter housing

## Part Numbers:

RB6-22X28X18-WB (4.3 sq. ft. floor)  
RB3-30X24X18-WB (5 sq. ft. floor - 3-Cage Rack)  
RB6-30X24X18-WB (5 sq. ft. floor - 6-Cage Rack)

## Dimensions:

61" W x 28.5" D x 75" H  
36"W x 30.5"D x 75"H  
64.5"W x 30.5"D x 75"H

## Designed for use with:



rabbit



scan here



more info

## Standard Features Include:

- Made from 304 stainless steel including the frame, feeders and pans
- 6 individual cages with 6 square feet per cage
- Autoclavable
- Side and back panels have solid lower with wire rod upper section for more light and visibility with better ventilation
- J feeders and stainless steel litter pans
- Removable polypropylene floors for animal comfort and quieter housing

## Part #

RB6-30X29X18-SS

## Designed for use with:



rabbit



## Standard Features Include:

- Made from PVC coated stainless steel mesh
- Two removable cages per unit
- Unit back made from plexiglass to protect facility walls
- Entire front panel is removable with smaller swinging style door that offers access inside the cage
- Two-handled crowding wall (with removable handles) to help capture the tree shrews

## Part #

TS2-20X24X80-C1

## Dimensions:

20" W x 24" D x 80" H

## Designed for use with:



tree shrew

# Glove Port Animal Isolator



## Standard Features Include:

- 304 stainless steel frame with large polycarbonate windows
- Includes both negative containment and positive isolation capability
- HEPA filtration for both supply and exhaust
- 3-speed fan control
- Sterile access port on one side
- Automated water system with adjustable height drinker

## Part #

AIU-01-PN-E

## Dimensions:

48" W x 28.5" D x 75" H

## Designed for use with:



rabbit



swine



## Standard Features Include:

- Made from all 304 stainless steel construction
- 5" stainless swivel casters with grease zerks, 2 with brakes
- Stainless steel excreta pans
- Side wall swings down for entry ramp
- Ramp has durable rubber matting and ridges for traction
- Grid wire top that opens
- Grid wire/solid stainless walls for greater visibility

## Part #

LGA-50.1X27.5X46.6-TRANSPORT

## Dimensions:

50.1" W x 27.5" D x 46.6" H

## Designed for use with:



canine



swine



sheep



cattle

# Open Top Animal Transport Unit



## Standard Features Include:

- Made from 304 stainless steel tubular, c-channel and rod construction
- 2 swing-out doors, one on each side of the cart
- Gate style slide door latches
- Removable doors
- Vinyl flooring coated #9 Tenderfoot®
- 8" high step-up
- Removable excreta pan from either side

## Part #

LGA-24X54X36-TC

## Dimensions:

24.5" W x 56" D x 36" H

## Designed for use with:



canine



swine



sheep



cattle



## Standard Features Include:

- Open top animal transport cart designed for sheep, pigs and dogs
- Made from 304 stainless steel tubular, c-channel and rod construction
- 2 swing-out doors, one on each end of the cart
- Includes a squeeze bar with sliding adjustable drop for reducing cage size
- Vinyl coated flooring
- Removable excreta pan from either side

## Part #

LGA-24X54X36-TRANSPORT-C3

## Dimensions:

24" W x 54" D x 36" H



## Designed for use with:



canine



swine



sheep



cattle

Modular Ventilated Caging

Large Animal Runs

Bulk Trucks

Avian Caging

## Proven Products, Creative Solutions.

For over 27 years, our company has built a reputation for quality products and friendly service. Facilities worldwide have trusted us to care for their animals and help further their research. With an industry-leading warranty, we'll continue to stand behind our products.

Our entire team will work tirelessly to craft a custom solution that meets your needs, fits your budget and exceeds your expectations. We want to partner with you on your next project. Contact us today!

Scan Here

## Alternative Design Manufacturing & Supply

Where you get quality products and service for all your small & large animal housing needs.

Whether it's an animal housing design to maximize your space usage, or material handling equipment to improve workflow efficiency, our commitment is to work closely with you to deliver the animal room solution that best meets your facility's space, needs and budget.

### Warranty



Alternative Design Manufacturing & Supply offers a 5-year warranty on all stainless steel fabricated products. These products are warranted to be free from material and workmanship defects during this period. The warranty period will begin from either the date put into service, installation (if applicable) or 6 months from date of shipment, whichever occurs first.

- Individually Ventilated Caging (IVC) systems including rack frames and fabricated components
- All other stainless steel fabricated caging systems
- Fabricated material handling, storage and processing equipment made of stainless steel



Alternative Design offers a 2-year warranty on all Polysulfone plastic caging products and electrical components. These products are warranted to be free from material and workmanship defects during this period. The warranty period will begin from the date of delivery. Our Polysulfone products will have a 2-year pro-rated warranty\*. Please note – Polycarbonate caging is not recommended for autoclaving and any warranty will be voided if autoclaved.

- All electrical components
- Polysulfone plastic caging products or components

\* Plastic caging materials that show physical, chemical or heat damage will not be covered under warranty.



Alternative Design offers a 1-year warranty on all other products. These products are warranted to be free from material and workmanship defects during this period. The warranty period will begin from the date of delivery unless otherwise stated.

- All other products manufactured or sold by Alternative Design Manufacturing & Supply

Should a failure occur within any of the the above listed time frames, Alternative Design will have the option to repair or replace the defective product(s), part(s) or component(s). This warranty is extended only to the original purchaser of the product and excludes damage due to abuse, misuse, as well as normal wear and tear. All non-warranty repairs will include a 90-day warranty from the date of delivery.

For more information on our industry-leading warranties, please contact Customer Service.



### Warranty Repair or Service

It is the goal of Alternative Design to provide a quality product with full service warranty support. Ask about Alternative Design's Performance Maintenance Agreement program. You may obtain warranty repair or service by contacting our Customer Service Department.

Alternative Design Manufacturing & Supply, Inc.  
3055 Cheri Whitlock Drive, P.O. Box 6330  
Siloam Springs, AR 72761  
USA

800-320-2459  
479-524-4343  
479-524-4125 (FAX)

info@altdesign.com  
www.altdesign.com  
f t in g+ You Tube